



## Contact us

+31(0)72 71 12 655 • [info@cellvision-technologies.com](mailto:info@cellvision-technologies.com)  
[www.cellvision-technologies.com](http://www.cellvision-technologies.com)

## Distributors worldwide



## Your local distributor

## Products



CellVision - Type Description	Articlenumber
-------------------------------	---------------

### 10 micron (µm) chamber height

CellVision, 2 chamber - horizontal	CV 1010-2ch
CellVision, 2 chamber - vertical	CV 1010-2cv
CellVision, 4 chamber - horizontal	CV 1010-4ch
CellVision, 4 chamber - vertical	CV 1010-4cv
CellVision, - 10x10 grid (Makler type)	CV 1010-102

### 12 micron (µm) chamber height

CellVision, 2 chamber - horizontal	CV 1012-2ch
CellVision, 2 chamber - vertical	CV 1012-2cv
CellVision, 4 chamber - horizontal	CV 1012-4ch
CellVision, 4 chamber - vertical	CV 1012-4cv

### 16 micron (µm) chamber height

CellVision, 2 chamber - horizontal	CV 1016-2ch
CellVision, 2 chamber - vertical	CV 1016-2cv
CellVision, 4 chamber - horizontal	CV 1016-4ch
CellVision, 4 chamber - vertical	CV 1016-4cv

### 20 micron (µm) chamber height

CellVision, 2 chamber - horizontal	CV 1020-2ch
CellVision, 2 chamber - vertical	CV 1020-2cv
CellVision, 4 chamber - horizontal	CV 1020-4ch
CellVision, 4 chamber - vertical	CV 1020-4cv
CellVision, 6 chamber - vertical	CV 1020-6cv
CellVision, 8 chamber - vertical	CV 1020-8cv
CellVision, - 10x10 grid (Makler type)	CV 1020-102
CellVision, - EU Post-Vasectomy grid	CV 1020-PV

### 30 micron (µm) chamber height

CellVision, 2 chamber - horizontal	CV 1030-2ch
CellVision, 4 chamber - horizontal	CV 1030-4ch

### 100 micron (µm) chamber height

CellVision, 2 chamber - horizontal	CV 1100-2ch
CellVision, 2 chamber - vertical	CV 1100-2cv
CellVision, 4 chamber - horizontal	CV 1100-4ch
CellVision, 4 chamber - vertical	CV 1100-4cv
CellVision, - Neubauer improved	CV 1100-NI
CellVision, - Bürker-Türk	CV 1100-BT
CellVision, - Thoma	CV 1100-T
CellVision, - Bürker	CV 1100-B

### 200 micron (µm) chamber height

CellVision, - Fuchs Rosenthal	CV 1200-FR
-------------------------------	------------

# A clear focus on counting cells



## Why Cellvision?


Easy to use  
Hygienic  
Accurate


applications

animal

industry

human

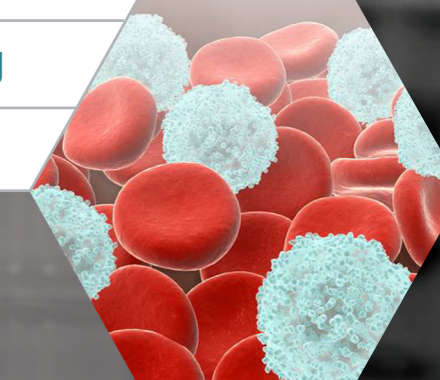

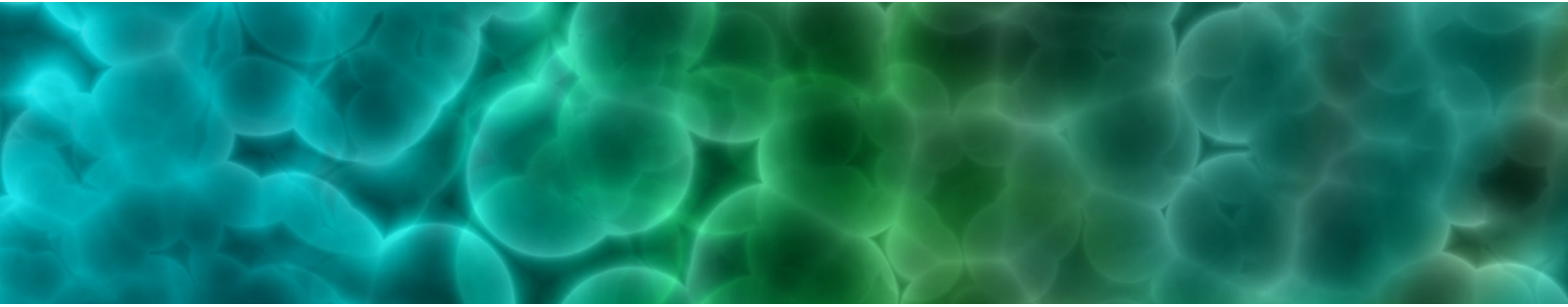

## Easy to use



Fill & Count glass slides with separate counting chambers and a fixed coverslip. CellVision slides have fitted counting grids with the exact same dimensions and identical image as re-usable hemocytometers. CellVision slides are also available without a built-in grid for automated cell analysis, cell culture and immunofluorescence assays. Designed for use with all cells/particles.

## Hygienic

Practical Infection Control: CellVision slides eliminate non-hygienic and time consuming cleaning procedures associated with re-usable hemocytometers - reducing the risk of exposure to potentially infectious material.


## Accurate

Precise fixed depth counting chambers with superior accuracy and reproducibility at a very comfortable price level.

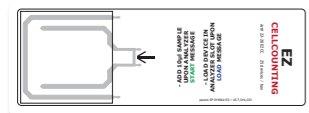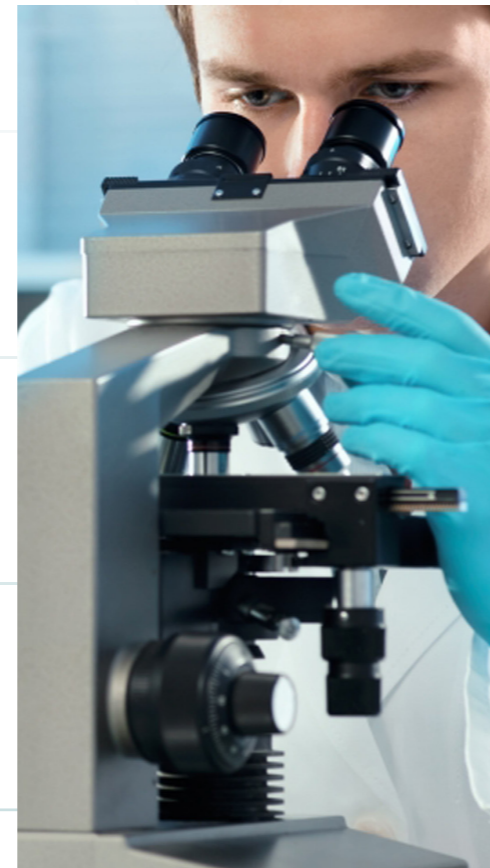

### EZ Cellcounting 22-20

Height	20 micron
Grid	n.a.
Chambers	1
Orientation	horizontal
Sample Volume	10 microliter
Shelf life	24 months

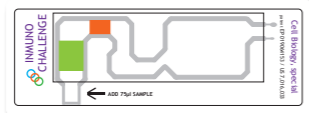

### Immuno Challenge - Cell Biology

Height	100 micron
Grid	n.a.
Chambers	1
Orientation	horizontal
Sample Volume	75 microliter
Shelf life	24 months

Custom designed slides are available upon request.

**Properties CV 1100-NI**


Height	100 micron
Grid	Neubauer Improved
Chambers	2
Orientation	horizontal
Sample Volume	25 microliter
Shelf life	24 months

**Properties CV 1100-B**

Height	100 micron
Grid	Bürker
Chambers	2
Orientation	horizontal
Sample Volume	25 microliter
Shelf life	24 months

**Properties CV 1100-BT**

Height	100 micron
Grid	Bürker-Türk
Chambers	2
Orientation	horizontal
Sample Volume	25 microliter
Shelf life	24 months

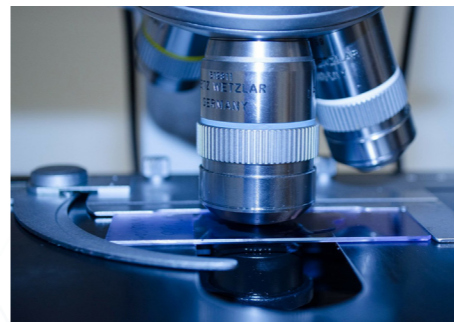

CellVision slides have many uses:

All common cell counts e.g. WBC, RBC, platelet, yeast, fungi, bacteria, spore, algae.

**Properties CV 1100-T**


Height	100 micron
Grid	Thoma
Chambers	2
Orientation	horizontal
Sample Volume	25 microliter
Shelf life	24 months

**Properties CV 1200-FR**

Height	200 micron
Grid	Fuchs Rosenthal
Chambers	2
Orientation	horizontal
Sample Volume	50 microliter
Shelf life	24 months

**Properties CV 1100-2CV**

Height	100 micron
Grid	n.a.
Chambers	2
Orientation	vertical
Sample Volume	25 microliter
Shelf life	24 months

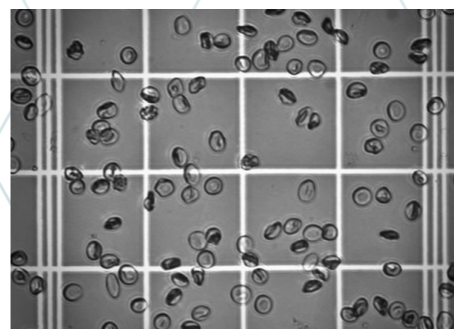

CellVision slides are designed to increase on accuracy, reproducibility, productivity and hygiene - saving on costs and time.